

CE

**14000**

**HIGH ELECTRICAL POTENTIAL  
THERAPY EQUIPMENT**

HPT003



**User Manual**



## Welcome to use our High Electrical Potential Therapy Equipment

Please read the manual carefully before operation, in order to guarantee the safe and full use of all the functions of this equipment.

# Index

I. Working Principle .....	1
II. Technical Parameter .....	1
III. Structure and Components .....	2
IV. Application and Contraindications .....	5
V. Contraindications .....	5
VI. Operating Instruction .....	6
VII. Energy Detection For Electric Field .....	12
VIII. Product Maintenance .....	13
IX. Product Service .....	18

# I. Working Principle

Electrical potential therapy equipment mainly adopts high-potential, negative potential, imitated natural electric field built by mid-frequency, to make the people under treatment back to the nature, repairing various of disorders and obstacles caused by the damages of natural environment.

After processing by this equipment, the alternating current (Ac) is transformed to AC high voltage alternating electric field, which can make all parts of the user vibrated subtly, therefore the imbalance of the body tissues and organs are fully recovered, the inherent healing capacity of human is improved, the cells will be activated effectively, the nerve conduction function will be strengthened, and the gastrointestinal motility will be accelerated so as to achieve the treatment of chronic constipation. By means of artificial method, the electric field with high vibration frequency is acted to the human body, spreading around the user and generating stable mid-frequency electromagnetic field of ac high-voltage negative electricity so as to promote blood circulation, unclog the collaterals, release pain and etc.

When using negative high voltage output, the negative high voltage electric field will break the original charge distribution of body, the changed charge swims among the tissues and the ceels then results in a series of biological or physical chemical changes which will promote the biological functions and pathological state of the tissues, and organs, makes parasympathetic excited to regulate vegetative nerve functional disturbance thereby treat neurasthenia, headache, insomnia ect.

# II. Technical Parameter

Power input:	AC 220V 50-60Hz
Input Power:	85VA
Fuse:	T1.6AL/250V
HPOT output voltage:	9000V 12000V 14000V
HPOT output frequency:	50Hz
Negative potential output voltage:	-6000V -8000V - 10000V
Negative potential output frequency:	50Hz
IFelectric potential output voltage:	2800Vp-p
IFelectric potential output frequency:	80 KHz
Display of reading :	LED Digital Screen
Default treatment time:	30min
Operating mode:	Continuous
Timing Range:	less than 90 minutes

Time Adjustment:	10min/gear, increasing circulation
Waterproof Level:	Non-watertight
Working environment:	Ambient temperature: 5°C ~ 40°C Relative humidity: ≤80 % Atmospheric pressure: 860-1060 hPa
Storage Environment:	Temperature range: - 40°C ~ 55°C Relative humidity: ≤80 % Atmospheric pressure: 500hPa~1060hPa
Host size:	430mmx270mmx331mm

# III. Structure and Components

## 1. Structure Diagram

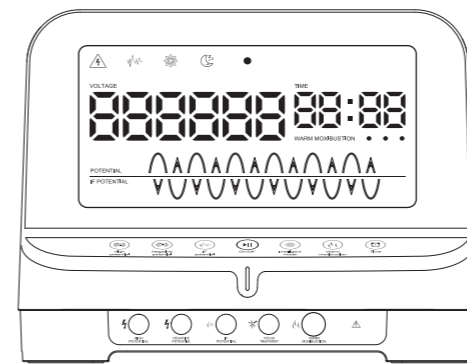


Figure 1

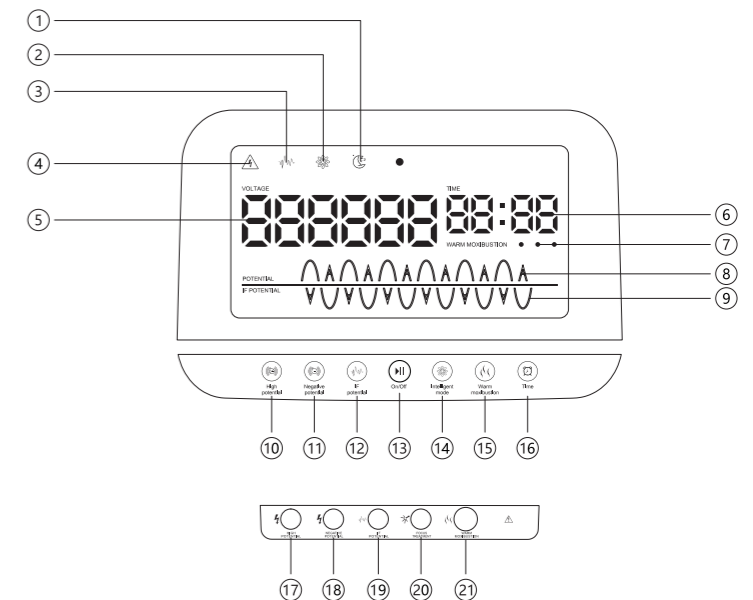


Figure 2

## 2. Structure Description (According to Figure 2)

1. Sleep mode icon display: When entering sleep mode, the corresponding sleep icon will light up.
2. Smart mode icon display: When entering smart mode, the corresponding smart mode icon will light up.
3. Intermediate frequency mode icon display: When entering the intermediate frequency treatment mode, the corresponding intermediate frequency icon will light up.
4. Voltage icon display: When entering high/negative potential treatment, the corresponding voltage icon will light up.

5. Voltage output display: Select different modes and gears, and the corresponding voltage values will be displayed; the default voltage at startup is 9000V.

6. Time display: The default display is 30 minutes when turned on. After starting treatment, the time will count down in seconds (sleep mode counts down in minutes); if the treatment is paused and switched to other modes or gears at this time, the time will accumulate and will not be reset. It is 30 minutes, but the time can be re-adjusted (adjustable from 10 to 90 minutes, in steps of 10 minutes).

7. Warm moxibustion mode display: When the warm moxibustion mode is selected, the corresponding icon will light up;

8. Intermediate frequency voltage output: During intermediate frequency treatment, the corresponding waveform will be displayed and "IF POTENTIAL" will light up at the same time.

9. High/negative potential voltage output: When entering high/negative potential treatment, the corresponding waveform will be displayed and "POTENTIAL" will light up at the same time.

10. High potential button (⊕): Press the first time to select the 9000V; press the second time to select the 12000V; press the third time to select the 14000V.

11. Negative potential button (⊖): Press the first time to select the -6000V ; press the second time to select the -8000V; press the third time to select the -10000V.

12. IF button (⊕): After pressing, the IF icon lights up, the corresponding voltage output waveform lights up, and "IF POTENTIAL" lights up at the same time.

13. On/Off key (⊕): start and stop button;

14. Smart mode button (⊕): Press it to turn on the smart treatment mode. The default is 9000V. After 5 minutes, it switches to 12000V. After another 5 minutes, it switches to 14000V. This cycle continues.

15. Warm moxibustion mode button (⊕): Press it to turn on the warm moxibustion treatment mode, which needs to be used with the corresponding warm moxibustion cushion.

16. Time button: adjust the time, adjustable from 10 to 90 minutes, in steps of 10 minutes.

17. High potential output port: Used in connection with the black energy pad, and the corresponding gear switching output is performed through the high potential button;

18. Negative potential output port: used in connection with the black energy pad, and the corresponding gear switching output is performed through the negative potential button;

19. Intermediate frequency output hole: used in connection with the black energy pad, and output by controlling the intermediate frequency button;

20. Pain releasing handle output jack: Used in connection with the Pain releasing handle, when using the high potential and negative potential functions, local treatment can be performed.

21. Warm moxibustion output port: used in connection with the warm moxibustion cushion, and output by controlling the warm moxibustion button.

### 3. Intelligent function output instruction

Function	Intelligence		
High potential circulation	9000V, 5 minutes	12000V, 5 minutes	14000V, 5 minutes

### 4. Sleep mode output description

Function	Intelligence			Default 1H, adjustable from 1-12 hours;
High potential circulation	9000V	12000V	14000V	
Negative potential circulation	-6000V	-8000V	-10000V	

#### Unsealed instrument

Taking out of the instrument and accessories from the carton, carefully checking whether there is any damage for each part during transportation and examining if items listed on the packing lists are complete. If there is any damaged or missing, please contact the dealer.

### 5. Standard Accessories

Name	Quantity
host	1 pc
Energy effect cushion	1 pc
Remote Control	1 pc
Pain releasing handle	1 pc
Power cord	1 pc
User manual	1 pc



Energy effect cushion



Pain Releasing Handle



## IV. Application

This product is suitable for the treatment of headaches, insomnia, bone and joint pain and chronic constipation.

## V. Contraindications

**Users who are using the following medical electrical equipment can not use this therapy equipment for treatment simultaneously.**

- Pacemakers or any other electronic transplant equipment.
- Artificial heart-lung machine which is a life-support system
- The ECG measuring equipment which taking along.

Taboo	<p><b>Users who are using the following medical electrical equipment can not use this therapy instrument for treatment simultaneously.</b></p> <ul style="list-style-type: none"> <li>★Pacemakers or any other electronic transplant equipment</li> <li>★Artificial heart-lung machine which is a life-support system</li> <li>★The ECG measuring equipment which taking along.</li> </ul> <p style="text-align: center;">— May cause serious problem —</p>
-------	---

**The following patients should be agreed by the professional doctors before using the equipment.**

- Acute disease patients
- Cardiopath
- Malignancy
- Patients with a high fever
- Pregnancy
- Being treated or people who are not feel well

**People within the metal into body are forbidden to use this product; During the onset of life-threatening disease should stop using this equipment.**

## VI. Operating Instruction


### Confirmation items before using

Must be under guidance by professionals when using the equipment.

- It must ensure that the power is 220V 50-60Hz.
- It must ensure that the host enclosure, power cord, energy pads and connecting lines are intact and no cracks, no damage.
- It must ensure that the power cord, energy pad and connecting lines are properly connected and the equipment can work normally.
- It must ensure that the host, the energy pads, connecting lines and the user are remain completely insulated state.
- Please do not use the product and other appliances on a same outlet.
- Please do not use the product under the environment of heat, humidity, dusty, smoky and vapor.



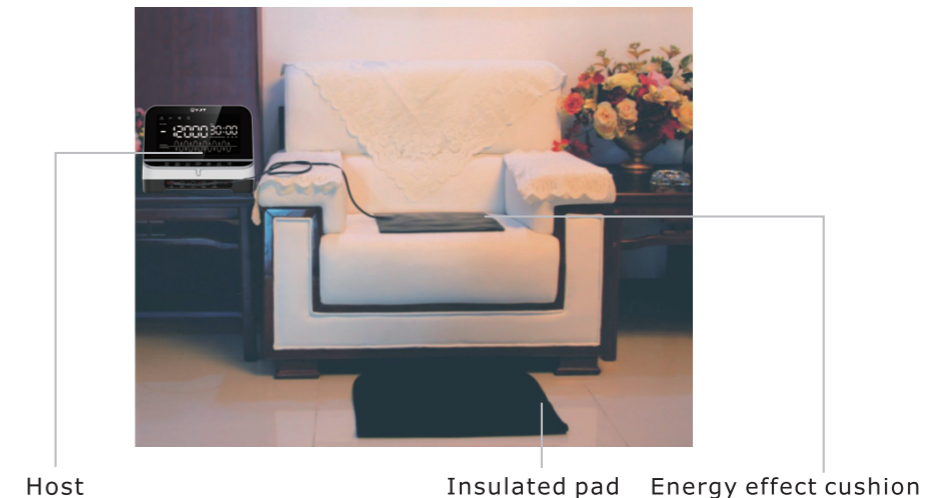

### Basic operation


Please operate in strict accordance with the following methods to ensure safe and proper use of the equipment.




- (1) Emplacing hosts. Flatly placed the host on a ventilation, dry easy-to-operation objects (plastic or wooden tables and chairs). The host must maintain a distance of at least 20cm with the ground the wall. Host and energy pad must maintain a distance of at least 1 m with computers and other appliances. It should be placed with a insulated pad on the ground for foot placing. See diagram.
- (2) Connection energy pad and power source. Making energy effect cushion (mattress) set smoothly on the plastic or the dry wooden chair (bed). The energy pad and wire cannot touch with the ground, the spacing maintain above 20cm. Inserting the plug into 220V 50-60Hz socket.
- (3) Turning on the host power button, the screen lightens the demonstration completely, the host enters the readiness for action, the timer demonstration is 30:00 minute.


(4) Function choice:

- A Before choosing high potential function, inserting the energy effect pad plug into host's high potential output hole, pressing the high potential button 10 , the high potential output may choose the 9000V, 12000V, 14000V of three kinds of intensities.
- B Before choosing negative potential function, inserting the energy pad plug into host's negative potential output hole, press the negative potential button 11 , the negative potential output may Choose -6000V, -8000V, -10000V of three kinds of intensities.
- C Before choosing intermediate potential function, inserting the energy pad and effect pad plug into host's intermediate potential output hole, press the intermediate potential button 12 , intermediate potential output specification are 2800V p-p, 80KHz.
- D When using the intelligent treatment, making the energy pad inserts into the host's high potential output hole, pressing the intelligent button 14 , electric potential output intensity is made by 5 minutes gap intelligence circulation of output high potential.
- E If you want to use the sleep mode, please insert the energy pad into the high potential or negative potential hole of the host, select the corresponding high/negative potential function, and select the appropriate gears. At this time, long press the time button to enter the sleep mode. The screen time will be displayed as "01:00", which means it can be used continuously for 1 hour; the time is adjustable, and the adjustable range is 1-12 hours.

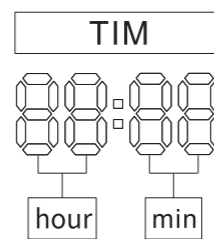

(5) Adjusting time.

Press the time key  to set the time, press to hand over by 10 minutes / grade for each time, from default 30:00, to 40:00, 50:00. ....Until 90:00 increasing. The most greatly fixed time for the count down time is 90 minutes. When count down, the time demonstrates for the first two digital presentations are minutes, the latter two digital presentation are seconds.

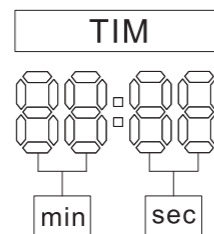

- (6) A. After setting treatment function and treatment time, after pressing start/stop key, the machine starts to work, the timer count down, after the period of using arrives, the buzzer sounds 8 times to remind the treatment has completed and withdraw from this treatment function.
- B. From step 4 to operate if the user would like to operate continuously.

(7) The Warm moxibustion function only can be use separately.

(8) After using, please switch off the host power, separate the power socket and make the host at the power-off condition.

## Items for operation reminding

- Laying a piece of insulated cloth when using the equipment on bed or sofa which attached metal spring.
- It must ensure that energy pad high pressure input line is inserted into high potential output hole when using high potential key.
- It must ensure that energy pad high pressure input line is inserted into negative potential output hole when using negative potential key.
- It must ensure that energy pad high pressure input line is inserted into intermediate potential output hole when using intermediate potential key.
- During use process, if transforming the function between the electric potential output and intermediate frequency electric potential outputs, may first press the key of On/Off to regain the readiness for action, then operating from elementary operation of 4th step to carry on the operation.
- During the use, if stopping in advance, it must press the key of On/Off to regain the readiness for action, then switch off the power.
- Users should drink 250-500ml lukewarm water before and after using the machine to enhance organism metabolism function and carry out waste.
- There are metabolic promotion responds during initial use, such as ache, dizziness, vertigo and a parched mouth. That is the normal phenomenon. Washing the face with the clear water after use, or the supplement of Vitamin ( eats fruit and vegetables, or oral administration of Vitamin C ). Symptoms will eliminate several days later.
- Heart disease, hypertension and diabetic's users: Starts to use one time for every day, starts from 30 minutes first, gradually increases to 40 minutes one month later. During use, please suspend for several days for later use again if has illness.

- For serious heart disease, hypertension and diabetic patients, it should be under the specialists instruction during the initial period use. The medicine cannot be taken decrement or stopped, if has ill please suspend using.

## Attention items during use

(1) In order to avoid the machine damage and an electric shock accident, please note:

- Do not place the machine to a place that the baby is easy to reach.
- During use, do not touch other people with the hand and other conductors of surrounding, do not touch the host's outer covering, do not transmit the thing mutually.
- The special attention to prevent the liquid to splash into the host and on the energy pad in order to avoid leaking age or damaging machine.
- Do not use the machine in the environment of moist, pollution and vibration.
- When treating in the winter and lying down on the bed, it cannot use with the electric blanket.
- When thunders, earthquake, fire, floods and power cut, it must immediately shut down the power source, captures the power line and suspend using.
- When the machine works, please do not hold the finger or put the metal goods and conductor insert into function output hole. If the high potential or negative potential function is not used, please plug the output port up with a plug.
- Do not serve as this machine outside the healthcare uses arbitrarily.

(2) For the following situation ,please use after the body returns to normal:

- After drinking wine excessively
- After strenuous exercise
- palpitation exceptionally
- scant of breath

(3) For the following situation, it should use under nursing of family member, please do not use seperately:

- Babies and infants
- Old and weak people
- Moving with difficulty
- Has epilepsy history
- Has the mentia and sensory disturbance
- Just start to use and does not understand application method
- Has aerosil and other filling materials in vivo
- Has the dysphrenia or the nerve medical history

## Time management

- It should be according to individual body adaptation to proceed in an orderly way and select the output function and service condition reasonably when starting to use the electric potential treatment meter. This product design default treatment time is 30 minutes, uses each time for 30 minutes. Only needing consistently insist that can make the effect. The physically weak person should reduce number and time of treatment, the strong person may relative extend time.
- When using intermediate frequency electric potential function, it should not exceed 60 minutes for each time. When the electric potential output is higher than 9000V, it should not exceed 60 minutes for each time. ( May refer to time management recommendation table )

The effect is better if twice for each day ( respectively one time for daytime and night). Uses for one month later. Adjusting the electric potential output intensity every 2-3 days. Looping use three grades output ; Or uses the intelligent key frequently , makes the circulation of upscale, middle -grade and low -grade. This may eliminate human body's tolerance to some constant electric field.

## Time management recommendation table

Function	1-10 days	11-20 days	21-30 days	30 days later
High potential	30min/ 9000V	30min/ 12000V	30min/ 14000V	Using the intelligent key, or makes the circulation of upscale, middle -grade and low -grade for every 2-3 days.
Negative potential	30min/ -6000V	30min/ -8000V	30min/ -10000V	
IF potential	30minutes for each time	30 minutes for each time	30 minutes for each time	Less than 60 minutes for each time

## Pain releasing handle' s function and operation tips

### Pain releasing handle' s function

Reducing nerve's excitability mainly through the high-voltage micro current stimulation. It can quickly reduce tingling, swelling, pain and other symptoms in some body parts.



### Potential strength and the use of time reference

Reasonably adjusting time and strength of potential as for the body's situation and recommending potential strength output is 9000V, 1-2 times daily, 10 minutes for each time. Diverting the water for the metal part of the dredge machine to some sensitive patients which feeling tingling so that to eliminate irritation.

### Operation tips


- Make the plug of pain releasing handle insert into host function output hole.
- Turning on the host power, setting using time.
- When the user is completely insulated, pressing the high potential/negative potential key and then pressing on/off button to enter the operation state.
- Holding the pain releasing handle insulation handles; aligned with the metal end to the affected area by pushing, clicking and massage.
- Neon tube will light up of the pain releasing handle's front-end when the pain releasing handle's metal side making effect on the human body. It shows that the machine is making effect on the patient's body part.
- If needing to stop in advance, first pressing on/off key to stop operation and recover to standby, then unplugging the machine. After using, please turn off the host's power and cut off the power switch to make the host in power off state.

### Notes for pain releasing handle's using

- Patients must use in a fully insulated state.
- Do not use the machine directly on the wound surface, ulcers and bleeding.
- The patient's hands, feet, body is not allowed to touch the wall, floor and the conductor during using the machine.
- Make sure the metal tip away from the computer or other electronic devices to avoid interference.
- Using the metal tip to make acupuncture point stimulation. It has a function of traditional Chinese medicine acupuncture.




### Diagrams of using methods

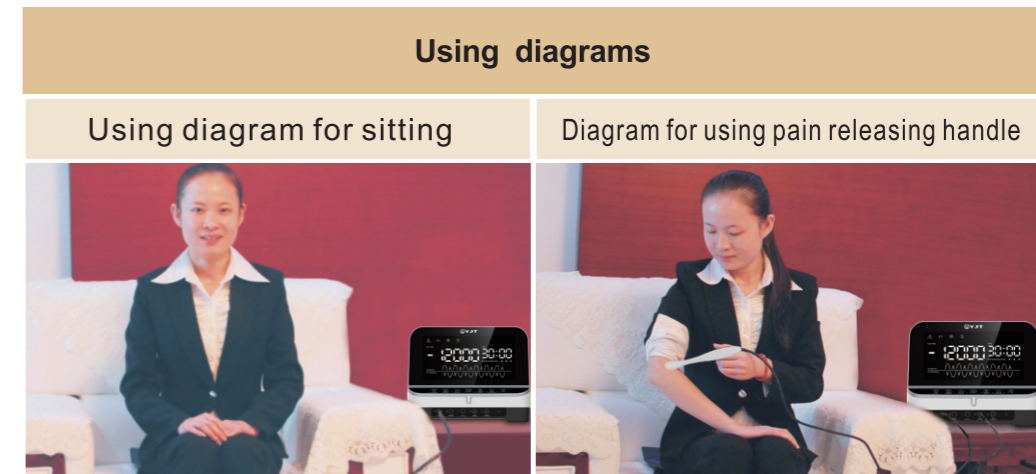

## VII. Energy Detection For Electric Field


Detector could check that if the host have an energy output for high potential, negative potential and IF potential. Also it could check if the user's body is in the range of above energy electric field. The detailed methods are as following:

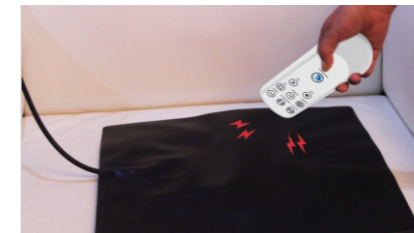

Holding the detector before using, making the detecting side near the sitting pad or lying pad which the energy is output. Pressing the detect switch, the detector's display lights up and beeps to show that the host and energy pad are in normal work state.



When using the machine, the user could hold the detector by oneself, make the detector side near the wall or other conductors except energy pad, then press the detector switch and the display lights up and beeps to show that the user's body is in the range of energy electric field.

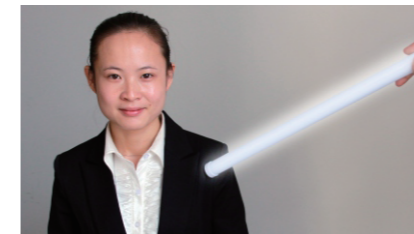

When using the machine, people outside the electric field could hold the detector, make the detector side near the user, then press the detector switch and the display lights up and beeps also to show that the user's body is in the range of energy electric field.

Please note: when using the machine, if the user holds the detector to oneself's body for checking, the detector display will not light up and not beep.

Note: the detector display will not light up and not beep for the following 2 situations:



1. The energy pad's plug is not inserted into host's energy output hole or the host is not start to output.
2. The detector inside battery depletion. It needs to open the battery cover to change for new batteries. (2 dry batteries of No.7)

#### **Unique detecting method for IF potential field energy:**

Beside for using detector to check IF potential field energy as above introduction, it could also check the IF potential field energy by household fluorescent tubes, the detailed methods are:

People who is outside the electric field holding the fluorescent tube by right hand, put the other side of the fluorescent tube on the user's body to keep touch with the user's skin, then stroking the bottom of fluorescent tube by left hand, the machine works in normal state if the fluorescent tube lights up. Colorful tubes can be used to facilitate the observation.

## **VIII . Product Maintenance**

### **Routine maintenance**

1. Do not place the instrument under the sun's direct rays, humidity, dusty and smoky environment;
2. Keep tidy for the appearance of the treatment instrument. It needs to be power off state when daily cleaning. It is not allowed to use organic solvents or chemical reagents to wipe the equipment. It can be wiped with a damp cloth. Make sure that the equipment is thoroughly dry before turning on and using;
3. When the product is not used, it should be kept in a cool, dry, ventilated and non-corrosive environment, and the power should be set up to check for once a month;
4. Move the host gently, it is not allowed to upside down, also not too hard when using keys, plugging and unplugging the connection line;
5. Make sure the energy pad and other accessories are good before using the machine. If there are defections, please contact the distributors or after-sale service for new replacement in time. Do not repair or disassemble the product by oneself, otherwise we can't provide a normal after-sales service for you.
6. Use soft cloth and concentration 75% medical alcohol to clean polluted energy pad is suggested. Do not use Benzene, gasoline and other chemicals to clean the machine and energy pad.
7. It should unplug the host line timely. Do not use wet hands to plug and unplug the power cord after using the machine for each time.
8. Non-folding and twisting high-voltage cable and power cord.

9. It is not allowed to touch the power switch and host switch by wet hands.
10. High-voltage cables should be firmly connected with the host interface.
11. Forbidden the foreign body enters to the host.
12. Immediately stop using for the following abnormal conditions:
  - a. Generation of smoke and unusual odor;
  - b. Water and other foreign matters entering into the host, high voltage cable and inside of the insulated cushion;
  - c. Generation of heat by the host, high voltage cable, insulation cushion and their joints;
  - d. Abnormal sound and lights of thunder, lightning and a lightning storm;
  - e. When the host, high voltage cable and insulated cushion are damaged.
13. It is required to turn off the power of the host before plugging high potential therapy instrument's accessories, moving the product and interrupting use.
14. Insulation cushion can not keep touch with the host after the host in the state of turn on.
15. Human body can finish use and leave the operation area after cutting off the power of host completely.

#### **Warning:**

As a number of high-precision parts are used in this machine, please do not disassemble this product, or try to repair it by oneself; immediately stop using if the instrument is in malfunction. It is not allowed to disassemble the chassis for repairing. Customers should contact the manufacturer or the dealer to repair broken machine promptly.

## Common trouble and remedy

Symptom	Checking and confirmation	Remedy
No display or no operation response	If the power cord is good connected	Please plug the power switch to keep good connection
	If the host is power on	Put the power turning on
	If the fuse ( tube ) is burned out	Checking and replacing the fuse of T 1.6AL /250V
Having the "buzzing" sound after power on	If the high-voltage wire keep touch with the ground and wall	Please set up high-voltage wire
Normally display but no output	If the energy pad is good connected	Re-connecting energy pad after turning the power off
	If the energy pad is connected to the proper output hole	Re-connecting energy pad after turning the power off

## 【 Tag Description 】



1. Make sure the warning sign on the treatment instrument is kept in a clear state, also to replace not clear or damaged sign.
2. Please familiarize with the manual to familiarize yourself with the operation of therapeutic equipment before using. Do not allow unqualified personnel to operate the treatment device.
3. Noting these warning and attention symbols attached on the treatment instrument or in this manual. They remind that the wrong operation may be dangerous to humans, or damaged to equipment.

Before installation of the treatment instrument into operation, please read the following safety precautions and warnings in details.



1. The energy jacks of this treatment instrument for high-potential, negative potential, IF potential are output for high voltage. Making sure to turn off the energy output when pulling and plugging to avoid electric shock.
2. Cutting off the power supply when overhauling, only qualified personnel can operate the treatment device. These personnel must be familiar with all of the warning symbols. Proper installation, operation and maintenance can ensure safety and maintain equipment properly.

## 【 Machine body plate 】

Note: This equipment is class II insulation against electric shock which is functional earth ground terminal and not protective ground.

## 【 Environmental statement 】

Earth is our home. It is everyone's responsibility to protect our environment. Please do not discard arbitrarily for the host and accessories of this therapy equipment after use!

## IX. Product Service

Dear Customer :

Thank you for purchasing our products, "customer satisfaction is the only standard" is our company's service concept. To protect your legal rights and eliminate the worry about using, our company made the following commitments and provide you with comprehensive and standard warranty services in turn.

### **(A) The warranty range:**

When purchasing this product, buyer is promised to have the warranty service for the host and related accessories. Please consult the local service provider for more details.

### **Entrust repair within warranty**

please refer to "Common trouble and remedy" for processing firstly.

if you still can not rule out failure, please kindly contact our service department.

please provide maintenance service guarantee.

### **(B) The following conditions are not included in warranty, our company will charge on some extent for repair.**

- 1) One can not provide the product warranty document and warranty card.
- 2) The damage caused by careless or incorrect use.
- 3) The machine failure is due to disassemble by oneself.
- 4) Due to damage caused by natural disasters.
- 5) The user who expiring the warranty period will be charged for materials fees properly.
- 6) Failure or damage due to accidental factors or man-made factors ( scratch, soaking, secretly replacing the configuration, etc.) Instrument outer case, brochures, packaging and other accessories.

### **Repairing that expiring the warranty period**

When surpassing warranty period, if it can also continuously use after repair, we may carry on the repair according to the customer's request and charge the essential repair fee.

**Inform to the user :**

Please fill in the product maintenance list in detail after purchasing this product so that we could put on file reference and provide post-sale service. The user should take responsibility if underserved due to absence of this procedure. Please understanding!

Warranty Card		For customer			
User Name		User Tel		Person In Charge	
Product Name		Product Model		Serial NO.	
Purchase Date		User Add			
Failures Description				Time	
Maintenance Record					
Failures Description				Time	
Maintenance Record					



Product Certification

The product has successfully passed factory quality inspection, testing contents are as follows:

- Appearance and controlling key
- Output voltage
- Short-circuit current
- Pressure Tester
- Leak test

Inspector \_\_\_\_\_

Inspection date \_\_\_\_\_