

成品尺寸:125*125mm
展开尺寸:125*250
对折页, 8P , 157g铜版纸 双面哑膜



NANO BELT



User Manual


CAUTION:

1. Do not place or store appliance where it can come into contact of any liquids.
2. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to persons.
3. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
4. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
5. Do not immerse unit in water. Keep liquids away from ventilation ports, buttons and charging port.
6. Do not operate continuously for more than one hour. After one hour use, allow unit to rest for 30 minutes before re-using.
7. Unplug the unit after charging or prior to use.
8. Do not remove screws or attempt to disassemble.
9. This is not a toy. For adult use only. Do not use if injured and consult your doctor before using this.


PRODUCT SPECIFICATIONS:

- Size of product itself: (L)1450*(W)240*(H)40.40mm
- Net weight of product itself: 1307g
- Single package size: 275*273* 120mm
- Battery (QTY & Type): 4 pc & 18650 2600mAh
- Certifications: UL

PRODUCT OVERVIEW:

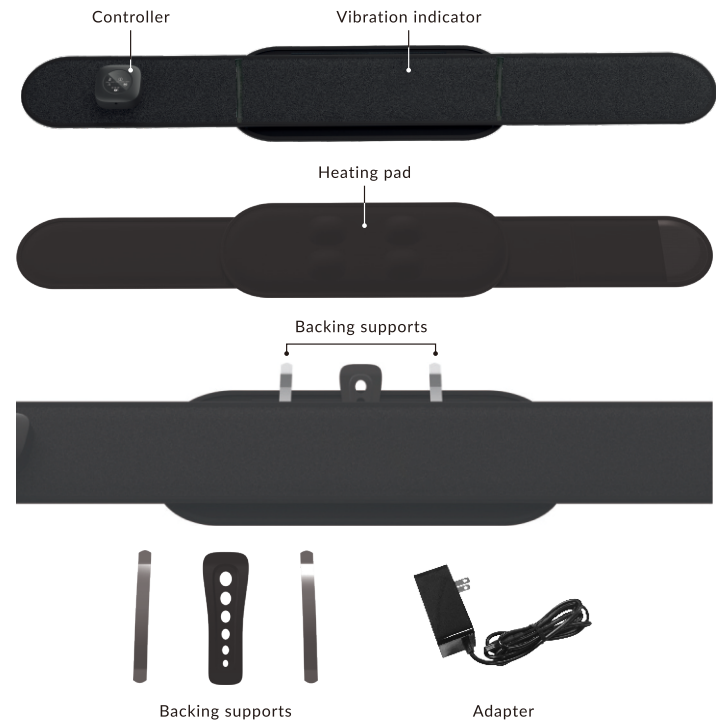

CONTROLLER INFORMATION:

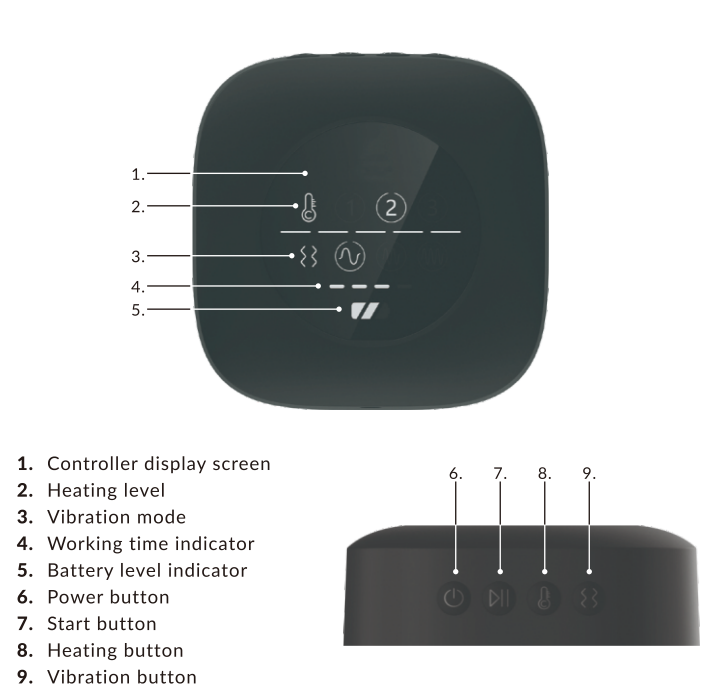

OPERATING INSTRUCTIONS:

- **Charging:**
 - To charge, insert the adaptor into the battery's charging port and plug into the charger into a wall outlet.
 - The LED lights of the battery level indicator (5) will flash to show the battery level and indicate active charging. Full charge is indicated when battery marks stop flashing.
 - Unit may be recharged at any time and at any battery level.
 - It is not recommended to fully deplete the battery to the lowest battery level. Recharge the unit monthly to best prolong the battery life.
- **Operating:**
 - Ensure battery is charged.
 - Wrap unit around the proper body part and attach straps.
 - Press and hold the power button (6) for 3 seconds to activate the unit. The screen will illuminate to display the default heating and vibration modes as seen in the image below.



- Press the start button (7) to activate unit.
- While using, press the heating button (8) to adjust heating level. By continuing to press the heating button (8) you will notice the LED screen heat level (2) will move between levels 1 through 3.
- While using, press the vibration button (9) to adjust vibration mode. There are three modes for your option. As you continue to press the button you will notice the three modes changing on the LED screen (3).
- Default working time is 20 minutes, and the unit will auto turn off after 20 minutes. The working time indicator (4) will indicate how far you are in the 20-minute treatment.

封底

封面

1

2

3

4

5

6