

# KNEE MASSAGER WITH HEATING



LED-K

User Manual

CE RoHS

CAUTION

- 1. Avoid looking directly at the light.
- 2. We do not recommend using this device if you have a pacemaker.
- 3. Please keep the light emitting ports clean to prevent light loss.
- 4. Improper use or adjustment of this device may invalidate your support services.
- 5. This device shall not be disassembled by yourself.
- 6. Please do not use this device under the environment of explosives or other inflammables.
- 7. You may feel hot while using this product, it's normal, please adjust the heating level when you feel too hot.
- 8. Please feel free to contact our customer service if there are any problems when using our product.
- 9. Keep children and pets away from the equipment. The equipment is designed for adult use only.
- 10. Always use the equipment as indicated. If you find any defective components whilst assembling or checking the equipment, or if you hear any unusual noise coming from the equipment during use, stop. Do not use the equipment until the problem has been rectified.


PRODUCT SPECIFICATIONS

- Battery (QTY & Type): 1 pc & Lithium battery 2450mAh
- Certifications: UL
- Size of product itself: (L)153.9\*(W)235\*(H)109mm
- Materials: ABS


FEATURES

- Physiotherapy
- Drug-free
- Chemical free
- Non-invasive
- Non-ablative
- Painless
- Requires zero downtime
- Safe for all skin types
- Safe for all ages
- No side effects
- Easy to use at home

INDICATIONS

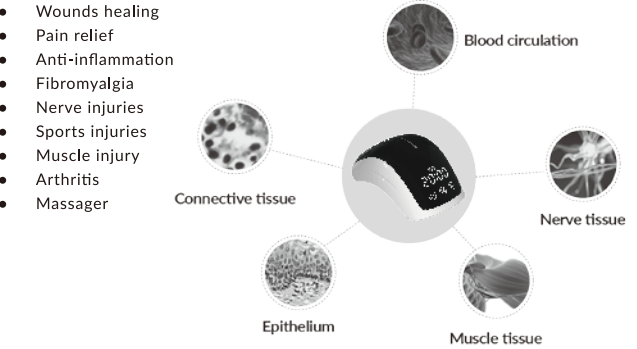

PRODUCT OVERVIEW

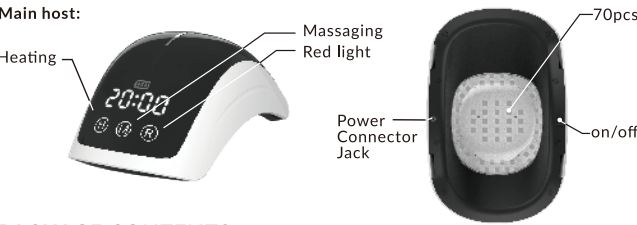


PACKAGE CONTENTS



OPERATING INSTRUCTIONS:

- Wear host on knee and tighten using strap.

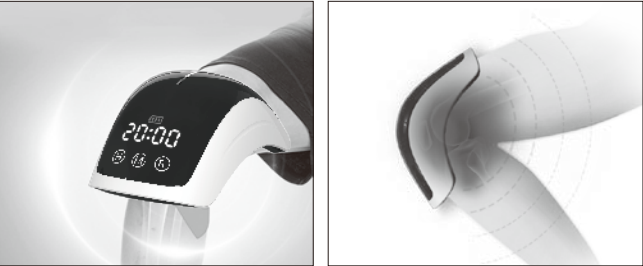

- **Turn On:**
  - Press on/off button at side of host device. The device will make a deep sound then it will turn on. The LED screen will display the default working time of 20 minutes.
- **Heating:**
  - Press H button to start heating function, there are 2 levels of heating intensity. The first time you press the H button, the power will go to lower intensity, when you press button again, the power intensity will begin the higher intensity function. To turn off the heating function press H a third time.
- **Massager:**
  - Press M button to start massaging function, there are 2 levels of massage intensity. The first time you press, the power will go to lower intensity, when you press button again, the power intensity will be higher. The vibrating function will turn off when you press M button a third time.
- **Red Light:**
  - Press R button to start red light and infrared light output, when you press button at second time the lights will turn off.
- **End of Treatment and Shutdown:**
  - All functions turn off after 20 minutes (default using time). You can also press on/off button to turn off the device to save battery.

TROUBLE SHOOTING

Defect	Causes	Solution
The button on the touch screen is not working	1. Battery power is not enough 2. Water got in, and the button is samaged	1. Please Charging immedately 2. The buttons are damaged that this device need to be returned for repairing
The light diodes is not working	Light diode is broken	Contact the agent

MAINTENANCE, CLEANING AND STORAGE

- 1. If a malfunction occurs, please do not try to disassemble this device yourself, if you cannot resolve the issue, contact an agent for help.
- 2. Please use cotton soaked in neutral cleaning fluid to clean the surface of this device. Do not use organic solvents, gas to clean it.
- 3. Please use cotton soaked in alcohol to gently clean light emitting ports. Do not use any liquid or hard objects to clean the light emitting ports.
- 4. Please keep this device away from humidity, high temperatures and direct exposure to sunlight to avoid damage or internal components aging before its expected lifetime.
- 5. Please do not place this device together with hard objects to avoid scratches on its surface or damage to the red lights.
- 6. General transportation tools are allowed, but strong shocks, vibrations, snow, or rain should be avoided in transportation process. Please store this device in a room without humidity or corrosive gases.